

CapellaScience

CapellaScience

precise

NEXT GENERATION
MICROPIPETTES



Experience truly
effortless
precision pipetting

The Capella Science range of Precise next generation micropipettes combines precision with ergonomics to deliver an effortless pipetting experience.



CapellaScience precise

NEXT GENERATION
MICROPIPETTES

Most current micropipettes are not up to the task of reducing fatigue or lessening strain on hands during use. Cutting-edge research requires equipment that is not only precise but efficient in design.

Reduced operating force ultimately leads to improved accuracy and repeatability. That's where the Capella Science Precise micropipette comes in. Its magnetic-assisted piston and ultra-low plunge force ensures your hands will not get tired even during extended bench work. Its advanced design provides maximum comfort and ease of use, thus delivering the perfect combination of precision and ergonomics.

Complete range of tips, from 10uL to 10mL.

Manufactured in state-of-the-art human-touch-free cleanroom, using medical-grade virgin polypropylene



Certified DNase, RNase, Pyrogen, Endotoxin, Human DNA & PCR Inhibitor Free.



Graduated



PE - filter tips for aerosol protection



Non-cytotoxic



Low-retention tips, free from softening agents



FDA-approved Medical grade PP




RACK PACK

- Fully Autoclavable & Stackable
- Hinge rack allows single hand operation
- Color coded tray for Easy identification





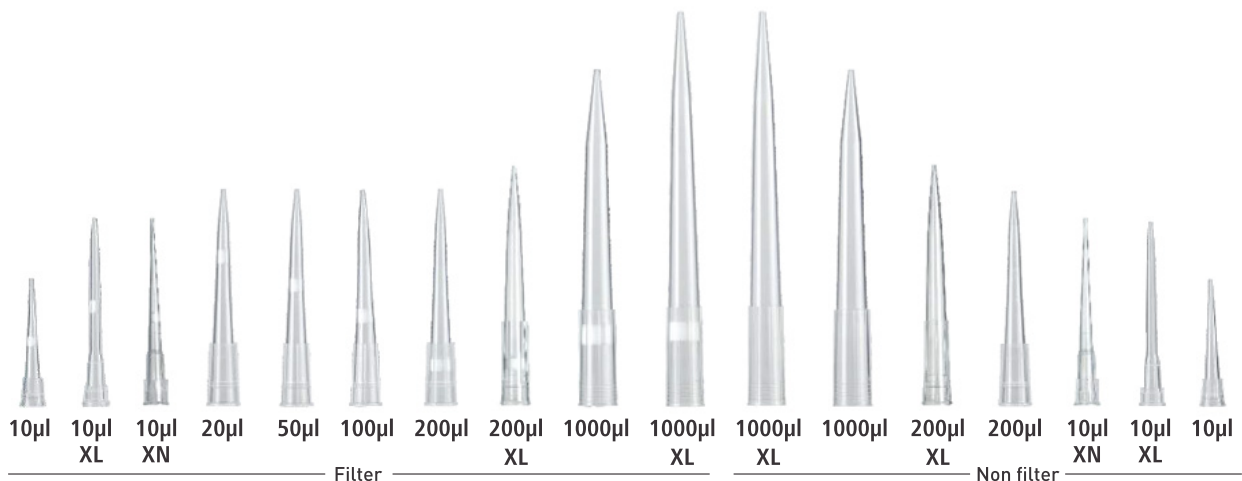


RELOADING STACK

- Reload rack with minimal effort
- Tips pre-arranged in 8X12 format
- Eco friendly with minimal plastic wastage




Full range of PCR Consumables, micro/centrifuge tubes, serological pipettes and deep-well plates also available, all manufactured in human-touch-free cleanroom conditions.

Extraordinary In Everything

-  Magnet assisted piston for consistently precise results
-  Large & clear 4-digit display
-  Made using highly durable material for protection against chemical corrosion & physical shocks
-  Spring-loaded tip cone for less strain & proper tip fitment for precise pipetting (multi channel only)
-  Highly durable PVDF tip cone
-  Color coding for easy identification
-  Ergonomic design with perfect weight distribution
-  Precise volume setting
-  3 year warranty
-  Innovative spring and seal design for minimal plunge force requirement
-  Calibrated as per ISO 8655-6 In ISO 17045 accredited lab
-  Easy in-house calibration & tool-free cleaning
-  Effortless Tip ejection with unique shock absorbing mechanism
-  Universal tip compatibility
-  Shelf mounting stand & carousel stand for easy storage
-  Durable Tip ejector made using special corrosion-resistant alloy
-  Volume lock setting
-  Fully autoclavable

Magnet assisted piston

Magnet assisted piston mechanism gives clear difference between first & second stop without using heavier spring. This assures consistent aspiration-dispensing with extremely low pipetting force

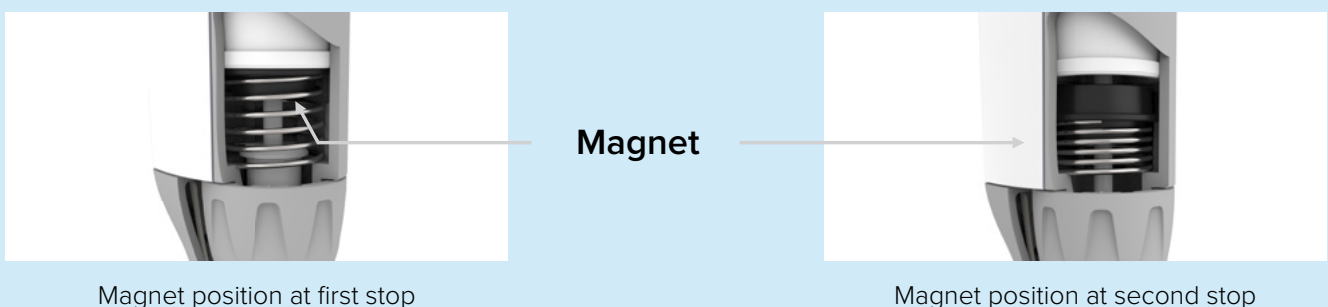

Ordering Information

CAPELLA SCIENCE PRECISE PIPETTES

Variable Volume			
Cat. No.	Volume Range [μ l]	Increments [μ l]	Introductory Pricing
NFV-2.5	0.1-2.5	0.002	\$155
NFV-10	0.5-10	0.02	\$155
NFV-20	2-20	0.02	\$155
NFV-50	5-50	0.1	\$155
NFV-100	10-100	0.1	\$155
NFV-200	20-200	0.2	\$155
NFV-1000	100-1000	1	\$155
NFV-5000	500-5000	10	\$165
NFV-10000	1000-10000	20	\$165

Fixed Volume		
Cat. No.	Volume Range [μ l]	Introductory Pricing
NFF-205	2.5	\$85
NFF-5	5	\$85
NFF-10	10	\$85
NFF-20	20	\$85
NFF-25	25	\$85
NFF-50	50	\$85
NFF-100	100	\$85
NFF-200	200	\$85
NFF-250	250	\$85
NFF-500	500	\$85
NFF-1000	1000	\$85
NFF-2000	2000	\$85
NFF-5000	5000	\$95
NFF-10000	10000	\$95

Multichannel: 8-channel			
Cat. No.	Volume Range [μ l]	Increments [μ l]	Introductory Pricing
NFM8-10	0.5-10	0.02	\$299
NFM8-20	2-20	0.02	\$299
NFM8-50	5-50	0.1	\$299
NFM8-100	10-100	0.1	\$299
NFM8-200	20-200	0.2	\$299
NFM8-300	30-300	0.2	\$299

Multichannel: 12-channel			
Cat. No.	Volume Range [μ l]	Increments [μ l]	Introductory Pricing
NFM12-10	0.5-10	0.02	\$345
NFM12-20	2-20	0.02	\$345
NFM12-50	5-50	0.1	\$345
NFM12-100	10-100	0.1	\$345
NFM12-200	20-200	0.2	\$345
NFM12-300	30-300	0.2	\$345

Pipette Stands		
Cat. No.	Volume Range [μ l]	Introductory Pricing
FS-01	Carousel stand for 6 pipettes	\$39
FS-02	Shelf mounting stand for Pipette	\$19



☎ 02 9575 7512

✉ enquiries@capellascience.com.au

🌐 www.capellascience.com.au